
ewimed

ewimed Pneu-Pack II

- DE** Gebrauchsanweisung
- EN** Instructions for use
- FR** Instructions d'utilisation
- IT** istruzioni d'uso
- SV** Bruksanvisning
- DA** Brugsanvisning
- NO** Bruksanvisning
- HU** Használati útmutató
- NL** Gebruiksaanwijzing

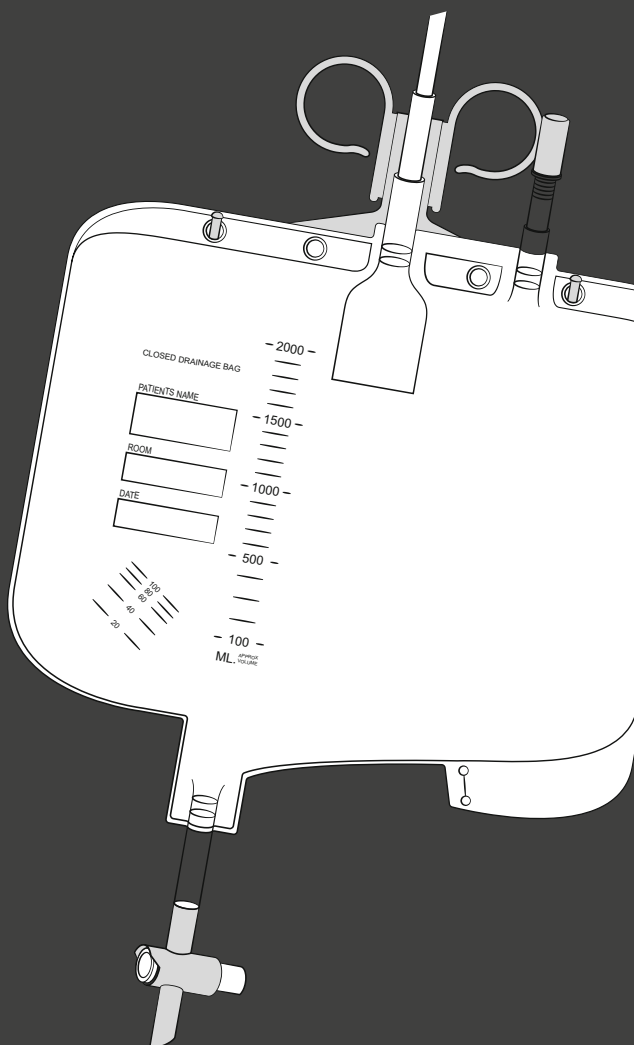REF P7899

DE ERKLÄRUNG DER SYMBOLE
EN EXPLANATION OF SYMBOLS
FR EXPLICATION DES SYMBOLES
IT SPIEGAZIONI DEI SIMBOLI
SV FÖRKLARING AV SYMBOLERNA

DA SYMBOLFORKLARING
NO FORKLARING AV -SYMBOLENE
HU SZIMBÓLUMOK MAGYARÁZATA
NL TOELICHTING VAN DE SYMBOLEN

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CONTENT

2	Explanation of symbols
8	Product description
9	Drainage Log Chart
22	The ewimed Pneu-Pack
22	Advantages
23	Warnings
	General
	In combination with the catheter
25	Connecting the Pneu-Pack II
27	Final steps
27	Glossary

EN

Read all the instructions for use carefully before using the product to drain an effusion.

All definitions of terms (*) can be found in the glossary on page 28.

YOU HAVE QUESTIONS OR NEED HELP?

Contact us:

DE: +49 7471 73972-0
info@ewimed.com

SE: +46 8 25 11 69
info.se@ewimed.com

BE: +49 7471 73972-0
info@ewimed.com

EN: +49 7471 73972-0
info@ewimed.com

DA: +45 31 43 50 60
info.dk@ewimed.com

NL: +31 203 695 620
info@ewimed.com

CH: +41 31 958 66 66
info.ch@ewimed.com

NO: +46 8 25 11 69
info.no@ewimed.com

LU: +49 7471 73972-0
info@ewimed.com

AT: +43 2231 22500
info@ewimed.com

HU: +36 1 998 8295
info@ewimed.hu

For further detailed information please visit our website www.ewimed.com

ewimed

THE EWIMED PNEU-PACK

The ewimed Pneu-Pack II is used for drainage and continuous drainage of effusions, such as pleural effusion* and ascites*, by gravity and for passive exhalation of air in pneumothorax* and pulmonary fistula*.

The Pneu-Pack II can be securely connected to the drainova®/drainova® ArgentiC catheters (art. no. 1140/1150/1160/1240/1250/1260) or the PleurX™/PeritX™ catheters (art. no. 50-7050/50-9050A) using the Luer-Lock female connector and a drainova® Irrigation Tube (art. no. 6060, 6060/2) or a PleurX™/PeritX™ rinsing tube/connection tube (art. no. 50-7245A).


The ewimed Pneu-Pack II is hereinafter referred to as the “reservoir”. The drainova®/drainova® ArgentiC catheter, PleurX™/PeritX™ catheters are hereinafter referred to as the “catheter”.

Minimum purchase quantity: 1PE PU (packaging unit) = 10 pieces, sterile

We also recommend the following for performing the drainage:
drainova® dressing set (art. no. 8020).


EN

ADVANTAGES

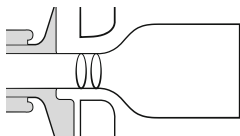
Integrated Heimlich valve* for maximum patient safety

The integrated Heimlich valve* ensures that effusions and/or air cannot flow back into the body once drained.

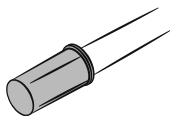
Air release valve

The air that collects in the reservoir, e.g. from a pneumothorax* or a pulmonary fistula *, can easily escape via the open air release valve.

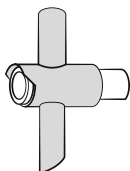T-tap valve

The effusion can be drained cleanly and regularly via the t-tap valve.

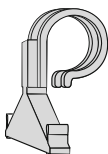Reservoir support

The reservoir support attached to the reservoir serves as a simple and secure hanging option.

WARNINGS

General

Do not use the product if:

- the sterile packaging is damaged
- the product is damaged
- the expiry date has passed

The reservoir is only suitable for use for the intended purpose.

The disposables are intended for single use only!

The rinsing tube is delivered in a sterile condition, which is assured by suitable quality controls. The reservoirs must not be resterilised.

The product is intended for single use only. Reuse may lead to cross-contamination.

ewimed GmbH assumes no liability for products that have been reused, resterilised or removed from damaged packaging.

Drainage from the chest cavity

Follow the doctor's instructions regarding the drainage quantity.

ewimed recommends draining no more than 1000 ml effusion per day. The drainage of larger effusion quantities should only be carried out on the instruction of the doctor treating the patient. The patient may feel an unpleasant pulling sensation during drainage. If this sensation becomes too intense, close the pinch clamp on the connection tube to pause the flow of liquid. Possible complications of recurrent drainage of a pleural effusion* include pneumothorax*, re-expansion pulmonary oedema*, hypotension*, circulatory collapse*, septa* in the chest cavity or infection*.

Drainage from the abdominal cavity


Follow the doctor's instructions regarding the drainage quantity.

ewimed recommends draining no more than 2000 ml effusion per day. The drainage of larger effusion quantities should only be carried out on the instruction of the doctor treating the patient. The patient may feel an unpleasant pulling sensation during drainage. If this sensation becomes too intense, close the pinch clamp on the connection tube to pause the flow of liquid. Possible complications of recurrent drainage of ascites* include hypotension*, circulatory collapse*, electrolyte imbalance*, protein loss*, hepatorenal syndrome*, leakage of ascites* from incisions*, peritonitis*, septa* in the peritoneal cavity* and infection*.

Position of the reservoir during drainage

As the effusion is drained by gravity, make sure that the reservoir is always situated below the catheter outlet during drainage.

Ensure that there is no tension at the catheter outlet and that the connection tube is not being pulled or tugged.

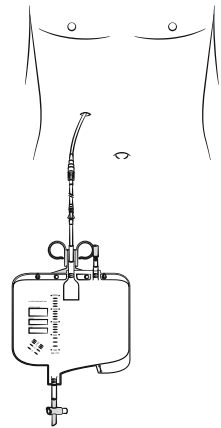Cleaning / disinfection

If effusion leaks and comes into contact with the skin, clean the skin with soap and water. For surfaces, use a suitable cleaning agent or disinfectant. Do not use octenisept® or disinfectants containing iodine to clean / disinfect the catheter or safety valve.

The alcohol swabs (wipes) are highly flammable.

Keep the swabs (wipes) away from naked flames.

Reservoir change intervals

The reservoir change should be carried out at a hospital or clinic following the usual hygiene standards. We recommend changing the reservoir every 72 hours at the latest.

The reservoir is CE-marked and sterilised with ethylene oxide gas (see packaging).

EN

WARNINGS


In combination with the catheter

Never cut through the catheter or cut off the safety valve.

Avoid direct contact between the catheter and sharp objects.

If the function of the safety valve is no longer guaranteed or if the safety valve has been accidentally cut off or is no longer present, proceed as follows:

- Firmly squeeze the catheter between your fingers.
- With your other hand, hold the emergency slide clamp open so that the catheter can be inserted into the open end of the clamp. It must be ensured that the emergency clamp is not applied directly to the valve, but instead at a distance of approx. 2 cm.
- Push the catheter fully into the narrow end of the clamp.

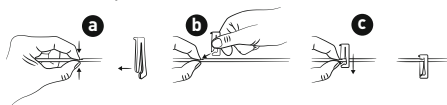**Only use original accessories.**

Only suitable original accessories may be connected to the catheter. Do not insert anything other than the access tip specifically designed for this purpose into the safety valve, otherwise the safety valve may be damaged. Improper handling can lead to air entering the body or to the continuous discharge of effusion.

Keep the safety valve on the catheter clean and the rinsing tube access tip on the connection tube sterile.


Make sure that these do not come into contact with non-sterile objects in order to avoid dirt or contamination. The safety valve must be disinfected before each drainage.

Make sure the safety valve and rinsing tube access tip are completely connected during the drainage process. Make sure the connection tube is not dragged or pulled.

When connecting the access tip to the safety valve, be sure to use a sterile procedure and ensure that the connection is firm. Accidental disconnection can result in contamination. In such cases, the safety valve must be disinfected with a fresh alcohol swab (wipe) and the disconnected drainage material must be disposed of. A new sterile reservoir must be used to continue the drainage. Please observe the local disposal regulations.

Note: If necessary, the plug-in connection can be secured with a separately available drainova® clickFix (art. no. P1100/P1100S) for additional protection against accidental disconnection.


The drainova® clickFix can be prescribed additionally and is not part of the drainage set.

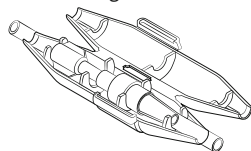
CONNECTING THE PNEU-PACK II

The Pneu-Pack II should be connected to the catheter following the usual hygiene standards at your institution.


1. Open the packaging, remove the reservoir and place it on a clean work surface. Make sure that the Luer-Lock connector on the connection tube remains sterile.
2. Remove a sterile drainova® Irrigation Tube or PleurX™/PeritX™ rinsing tube/connection tube from the packaging.
 - **drainova® Irrigation Tube:** Fully close the crimp clamp (A) on the connection tube.

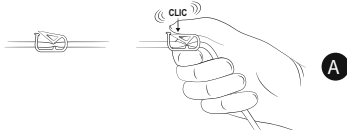- **PleurX™/PeritX™ rinsing tube/connection tube:** close the roller clamp (B) by turning it towards the Luer-Lock male connector.

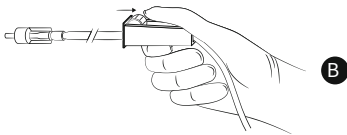
3. Connect the reservoir's Luer-Lock female connector to the Luer-Lock of the rinsing tube/connection tube (art. no. 6060,6060/2 (A), art. no. 50-7245A (B)).

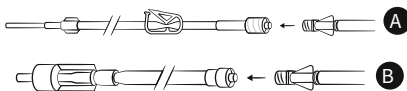
4. Loosen the drainova® silicone cap, sterile or the PleurX™/PeritX™ valvecap:
 - **drainova® silicone cap, sterile:** Firmly hold the catheter's safety valve and then carefully pull off the drainova® silicone cap, sterile. Dispose of the cap.

- **PleurX™/PeritX™ valve cap:** Firmly hold the catheter's safety valve and then twist the PleurX™/PeritX™ valve cap anticlockwise and gently pull it off. Dispose of the cap.

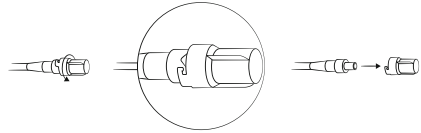
5. Hold the safety valve firmly in one hand and clean the area around the safety valve opening, e.g. with an alcohol swab (wipe) or a suitable alcohol-based disinfectant.

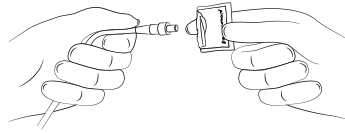Caution: Do not use octenisept® or disinfectants containing iodine.


- **drainova® Irrigation Tube:**

Hold the neck of the access tip and carefully pull the protective cover off the access tip, with a slight twisting movement if necessary.

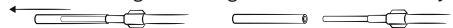Push the access tip into the safety valve.

Make sure that the safety valve and the access tip are firmly connected during the drainage process. You should hear and feel them snap into place.

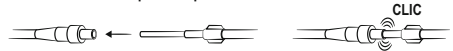- **PleurX™/PeritX™ rinsing tube/connection tube:** Firmly hold the cap and carefully pull the protective cover off the access tip.


Then push the access tip of the rinsing tube/connection tube into the safety valve and secure the connection by turning it. You should hear and feel them snap into place.

26 | Connecting the Pneu-Pack II

EN

6. Open the pinch clamp (A) / roller clamp (B) of the rinsing tube / connection tube to allow drainage to begin. Always keep the reservoir below the catheter outlet, as drainage is by gravity. The effusion now flows into the reservoir and/or can be used for passive exhalation of air in pneumothorax* and pulmonary fistula*. When draining ascites*, the roller clamp should be opened or closed at intervals according to the doctor's instructions.

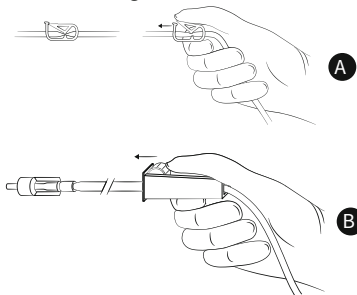Note: The flow of liquid may slow down or be interrupted if:


- the effusion is almost completely drained
- the effusion has a high viscosity
- the safety valve is blocked by a clot* or a protein thread.

7. The air release valve can be opened if required. This valve is closed on delivery.


Caution: If the reservoir is in a horizontal position, effusion may escape through the open air release valve.


8. To hang the reservoir, fold the reservoir support backwards and hang the support on, for example, the bed frame.

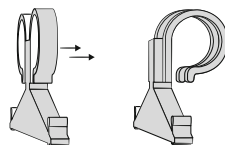9. This continuous drainage system is a safe and closed system. The reservoir and the rinsing tube / connection tube should be changed after 72 hours at the latest or in accordance with your institution's hygiene regulations.
10. To empty the reservoir, open the t-tap valve.


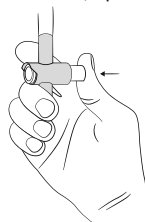
11. When the desired drainage volume has been reached, stop the drainage by closing the pinch clamp (A) / roller clamp (B).

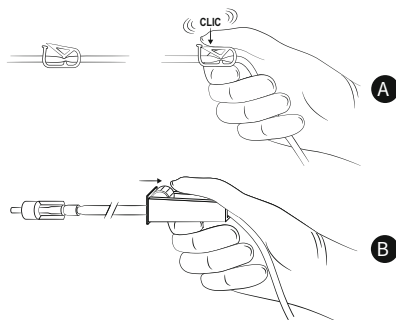FINAL STEPS


1. Open the packaging of the drainova® silicone cap, sterile or the PleurX™/PeritX™ valvecap but leave the cap in the packaging for the time being.
 - **drainova® Irrigation Tube:** Take the connection tube in one hand and the safety valve in the other hand. With a firm, steady motion, pull the access tip out of the safety valve. Remove the connection tube.
2. Clean the safety valve with an alcohol swab (wipe) or a suitable alcohol-based disinfectant.
3. Attach the drainova® silicone cap, sterile or the PleurX™/PeritX™ valvecap:
 - **drainova® silicone cap, sterile:** Push the new sterile drainova® silicone cap firmly onto the catheter's safety valve. The cap stays in position by itself without additional locking.

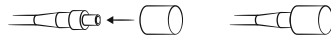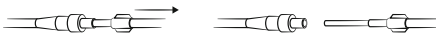- **PleurX™/PeritX™ rinsing tube/connection tube:** Take the rinsing tube/connection tube in one hand and the safety valve in the other hand. Loosen the connection with a twisting motion and pull the access tip of the rinsing tube/connection tube out of the safety valve with a firm, smooth motion.

- **PleurX™/PeritX™ valvecap:** Place the new sterile valve cap over the safety valve and turn it clockwise until it clicks into place.

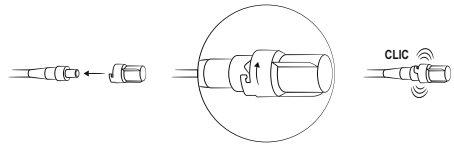EN

GLOSSARY

Ascites	Effusion or accumulation of fluid in the abdomen (stomach)
Circulatory collapse	Temporary loss of consciousness
DEHP	Plasticisers in plastic
Electrolyte imbalance	Disturbance of the balance of ions (e.g. sodium, potassium, calcium) in the blood or body cells
Heimlich valve	Protection against backflow for increased patient safety
Hepatorenal syndrome	Decrease in kidney function
Hypotension	Low blood pressure
Incision	Cut into the skin made during surgery, e.g. for implantation of the catheter
Infection	Reaction of the body's own defences to pathogens such as bacteria, viruses, etc.
Peritoneal cavity	Part of the abdominal cavity

Peritonitis	Inflammation of the peritoneum
Phthalate	Plasticiser (e.g. DEHP)
Phthalate exposure	Coming into contact with phthalates
Pleural effusion	Effusion or accumulation of fluid in the pleural cavity
Pneumothorax	Air in the space between the lungs and the chest wall
Protein loss	Reduction in the amount of protein in the blood (protein depletion)
Pulmonary fistula	Open connection between lung tissue and pleural cavity
Re-expansion pulmonary oedema	Retention of fluid in the functional lung tissue due to the sudden expansion of the previously compressed lung; acute dyspnoea
Septa	Tissue bridges/chambering between separated spaces in a body cavity (e.g. abdomen)
Thoracic drainage	Used to transport blood, secretions or air away from the pleural cavity (space between the pleura and the lung lining)

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PleurX™/PeritX™ is a trademark and/or registered trademark of Becton Dickinson and Company, or one of its subsidiaries.
drainova® is a trademark and/or registered trademark of ewimed GmbH.
octenisept® is a trademark of Schülke & Mayr GmbH.
ewimed Netherlands B.V.